



Clearly different.  
**Clearly better.**



Women's Health | HIV / STDs | Respiratory | Cardiovascular | Enterics & GI | General Health

## Clearview® hCG

### Clearview hCG Testing

One-step pregnancy  
tests for early diagnosis.

### Clearly different

*Clearview's patented innovative technology provides highly sensitive and specific tests for pregnancy.*

- **Rapid Response:** Single-step; positive results in as little as 1 minute.
- **Simplicity:** CLIA waived (urine), easy-to-use test.
- **Reliability:** Highly sensitive assay with built-in controls ensures accuracy.
- **Flexibility:** Use with multiple sample types — urine or serum.

### Clearly better

*Clearview provides better patient care with early diagnosis for proper planning, reducing potential risks to mother and fetus.*

- Earlier detection allows for improved care.
- Rapid results help prevent infection or contraindicated treatment.



**Clearview hCG Combo II**

40 Cassettes

Order No. 92215

Urine - CLIA Waived

Serum - Non-waived

**Clearview 25 hCG Combo**

40 Cassettes

Order No. 92225

Urine - CLIA Waived

Serum - Non-waived

**Clearview hCG Urine**

50 Dipsticks

Order No. 92211

CLIA Waived

**Clearview hCG Urine**

40 Cassettes

Order No. 92217

CLIA Waived

**hCG Controls**

1 Positive, 1 Negative

Order No. 92218

visit

[www.invernessmedicalpd.com](http://www.invernessmedicalpd.com)

for more information

**CLIAwaived.com**

For more information please contact CLIAwaived Inc. at  
 11578 Sorrento Valley Road ■ San Diego, CA 92121  
 Tel: (858) 481-5031 ■ Toll Free: (888) 882-7739  
 E-mail: [info@cliawaived.com](mailto:info@cliawaived.com)

Visit us on the web at: [www.cliawaived.com](http://www.cliawaived.com)  
 Distributed by: CLIAwaived, Inc.

**Clearview® hCG Combo Tests Comparison Grid**

Test Name	Clearview hCG Combo II & 25 hCG Combo	Quidel QuickVue® One-Step hCG Combo	Beckman Coulter ICON® 25 hCG	Genzyme OSOM® hCG Combo Test
Format	40 Cassettes	50 Cassettes	25 Cassettes	25 Cassettes
CLIA Status: Sample Types	Waived: Urine Non-waived: Serum	Waived: Urine Non-waived: Serum	Waived: Urine Non-waived: Serum	Waived: Urine Non-waived: Serum
Sensitivity	<b>hCG Combo II</b> 20 mIU/mL Urine 10 mIU/mL Serum  <b>25 hCG Combo</b> 25 mIU/mL	25 mIU/mL	25 mIU/mL	20 mIU/mL Urine 10 mIU/mL Serum
Specificity	>99%	>99%	>99%	>99%
Test Steps	1	1	1	1
Minutes to Result	(+) = 1* (-) Urine = 3 (-) Serum = 5	Urine = 3 Serum = 5	Urine = 3 Serum = 5	Urine = 3 Serum = 5
Shelf Life	24 months	18 months	24 months	24 months
Storage	Room Temperature	Room Temperature	Room Temperature	Room Temperature

**Clearview hCG Urine Dipstick Procedure**Dip to line,  
5 secondsRead in  
3 minutes


Interpret results



Positive




Negative

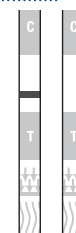

Invalid

**Clearview hCG Combo II / 25 Combo / Urine Cassette Procedure**

Draw sample

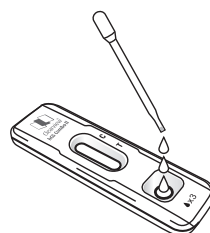Add 3 drops  
to sample wellRead in  
3 minutes - Urine  
5 minutes - Serum

Positive



Negative



Invalid

Interpret results

**CLIAwaived.com**

For more information please contact CLIAwaived Inc. at  
 11578 Sorrento Valley Road ■ San Diego, CA 92121  
 Tel: (858) 481-5031 ■ Toll Free: (888) 882-7739

E-mail: [info@cliawaived.com](mailto:info@cliawaived.com)Visit us on the web at: [www.cliawaived.com](http://www.cliawaived.com)

Distributed by: CLIAwaived, Inc.