

## Introducing our Exclusive Detector Portal



**Whidbey Island Clinic - Oak Harbor**  
 32550 State Route 20  
 Oak Harbor, WA 98277  
 Phone: (360) 555-5555 FAX: (360) 555-7777  
[contact@whidbeyclinic.com](mailto:contact@whidbeyclinic.com)  
<http://www.WhidbeyClinic.com>

Result	Patient Name	Date of Birth
COVID-19 Virus Not Detected	Anthony Simmons	01/20/1986

Sample Collected	Type of Test
2023 February 9, 20:55	SARS-CoV-2 (COVID-19) RNA [Presence] in Respiratory specimen by NAA with probe detection

Patient Identification		
Drivers License	Florida	G789-987-99-876-0

Test Device: The AscencioDx® Molecular Detector  
 Test Identifier: The AscencioDx® COVID-19 Test  
 Test Ordered: 94500-6  
 Test Collection Method: Anterior nares swab (697989009)  
 Test Collection Date: 2023 February 9, 20:55

Clinician: Linsy Collins, NP

Facility: Whidbey Island Clinic - Oak Harbor  
 Facility Address: 32440 State Route 20, Oak Harbor, WA 98277  
 Facility Phone Number: (360) 555-5555  
 Facility FAX: (360) 555-7777  
 Facility CLIA: 8798098079

Paula Hamshaw, MD  
 General Manager

Report letter generation date: 02/09/2023 21:05:16

- Comprehensive report generation and customizable "Passport" style report templates

PH

Start New Test
Scan Test Result

**Settings** | Organization Details | Patient Portal | Report Letters

Report Letters provide your patients with a formal copy of their test results for use as proof of negative COVID-19 test.

**Report Letter Preview**

**Options**

Template

DevTemplate.html (view sample or download)

DevTemplateGreen.html (view sample or download)

DevTemplateRed.html (view sample or download)

Upload New Template

**Options**

**Official Signature**

Letters must include the signature of an official representative of your organization. If you are the signee for your organization, please enter your name and title here to generate your signature. This name must match the name of Paula Hamshaw associated with this user account.

First Name:  Last Name:

Credentials:  Title:

*Paula Hamshaw*  
 Paula Hamshaw, MD  
 General Manager

Upload



Lit. Re-order No: MKTG-03  
 2023-03-23-V3.1.0

[Login](#) | [Create an Account](#)

The AscencioDx® Molecular Detector is a rapid, highly accurate, easy-to-use and affordable clinical diagnostic platform for performing molecular analysis of patient samples in point-of-care (POC) settings using The AscencioDx® COVID-19 Test.

The AscencioDx® COVID-19 Test and The AscencioDx® Molecular Detector has not been FDA cleared or approved, but has been authorized for emergency use by FDA under an EUA for use by authorized laboratories; This product has been authorized only for the detection of nucleic acid from SARS-CoV-2, not for any other viruses or pathogens; and the emergency use of this product is only authorized for the duration of the declaration that circumstances exist justifying the authorization of emergency use of in vitro diagnostics for detection and/or diagnosis of COVID-19 under Section 564(b)(1) of the Federal Food, Drug, and Cosmetic Act, 21 U.S.C. § 360bbb-3(b)(1), unless the declaration is terminated or authorization is revoked sooner.

[Learn more >>](#)
[Anavasi Diagnostics >>](#)

**Begin using your AscencioDx**


An Anavasi Diagnostics account allows you to manage patient records as well as documentation for operators who use your detector, review the testing history of patients, and simplifies the day-to-day use of your detector.

Login
Create an Account

Anavasi fully encrypts your data and securely sends it to the relevant federal, state, and local health authorities.

- Results are securely transmitted from the detector to our cloud-based portal via a QR-Code
- No Bluetooth required or app to download
- HIPAA compliant
- Automatically disseminates results to the CDC ReportStream, which in turn reports to in-program State and Local Health Agencies
- Intuitive and versatile dashboard that uses middleware to fully integrate with your existing LIMS and EMR systems
- Eliminate the hassle of paper or cumbersome spreadsheets while easily and securely managing patient records on your dashboard

- Set up your organization and manage detectors, locations, operators, physicians and patient data on your dashboard

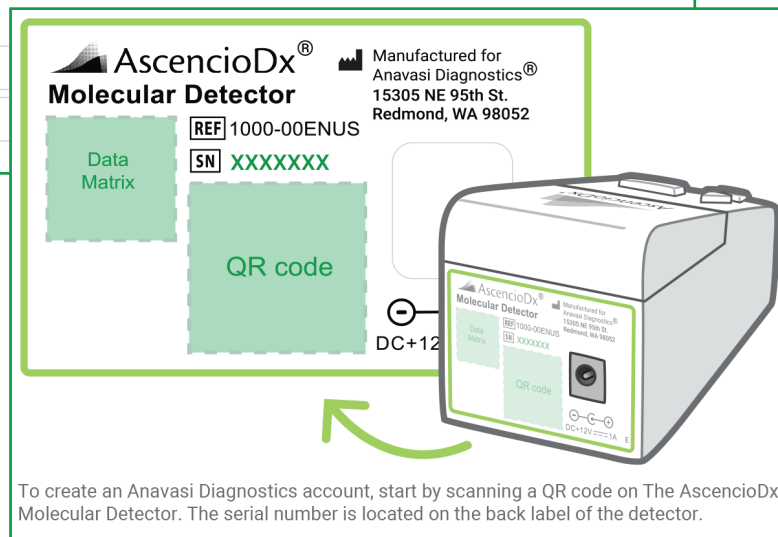

To create an Anavasi Diagnostics account, start by scanning a QR code on The AscencioDx Molecular Detector. The serial number is located on the back label of the detector.

- Value added portal exclusively available for The AscencioDx® Molecular Detector at no additional cost

PATIENT	TEST RESULT	DATE / DETECTOR	PATIENT PORTAL
<a href="#">Assign Patient</a>	POSITIVE COVID-19 Virus Detected	02/09/2023 12:30 pm [SN] 7385903 D1	<a href="#">Generate Letter</a> <a href="#">Email / Invite</a> <a href="#">Print / Download</a>
<a href="#">Assign Patient</a>	POSITIVE COVID-19 Virus Detected	02/09/2023 11:10 am [SN] 7385903 D1	<a href="#">Generate Letter</a> <a href="#">Email / Invite</a> <a href="#">Print / Download</a>
✓ <b>Jenkinson, Anthony Martin</b> DOB: 01 / 01 / 2001 @ ( ) -	POSITIVE COVID-19 Virus Detected	02/09/2023 6:40 am [SN] 7385903 D1	<a href="#">Regenerate Letter</a> <a href="#">Email / Invite</a> <a href="#">Print / Download</a>
✓ <b>Smith, Becky</b> DOB: 01 / 02 / 2001 @ ( ) -	NEGATIVE COVID-19 Virus Not Detected	02/09/2023 6:04 am [SN] 7385903 D1	<a href="#">Regenerate Letter</a> <a href="#">Email / Invite</a> <a href="#">Print / Download</a>
✓ <b>Skinner, Zak</b> DOB: 02 / 02 / 2002 @ ( ) -	NEGATIVE COVID-19 Virus Not Detected	02/09/2023 6:29 am [SN] 7385903 D1	<a href="#">Regenerate Letter</a> <a href="#">Email / Invite</a> <a href="#">Print / Download</a>
✓ <b>Houston, Daiana Amelle</b> DOB: 02 / 02 / 2002 @	POSITIVE COVID-19 Virus Detected	02/09/2023 6:09 am [SN] 7385903 D1	<a href="#">Regenerate Letter</a> <a href="#">Email / Invite</a> <a href="#">Print / Download</a>

- A click of a button securely emails test results to a patient
- Provide individual portal access to patients so they have a secure record of all their COVID results
- Results are also exportable in an excel format
- Patient information can be auto-populated by scanning their Driver's License – reducing data entry errors and saving time.